

Ease of use

No additional equipment needed for placement

One size fits all

For all gastric tubes up to 18 Fr

Economic

Quick placement of gastric tube saves OR time



Gastric Tube Guide

Insertion aid for the oral placement of gastric tubes

Gastric Tube Guide

Insertion aid for the oral placement of gastric tubes.
In order to avoid nasal injuries, the majority of gastric tubes that are not required postoperatively will generally be placed orally. For example, the gastric tube is inserted through the mouth after endotracheal intubation in patients who require a laparotomy or laparoscopy intervention without opening the intestine or stomach.


The introduction of a gastric tube is often more difficult and time-consuming than other anesthetic procedures (intubation, IV access, etc.).

The Gastric Tube Guide is the perfect tool for quick oral placement of gastric tubes. It is effective without complications.



Fast lane

Insertion funnel on proximal tip to ease the placement of all types of gastric tubes (PVC, polyurethane, silicone).

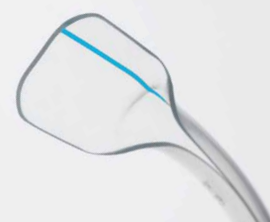

Design


Precurved shape to overcome the anatomical structure of the pharynx without the use of additional equipment.
Graduation marks in cm to adjust placement of the device.



Slit shaft to allow easy removal of the guide after placement of the gastric tube.

Atraumatic

Smooth and beveled distal tip to reduce risk of mucosal lesion during insertion.

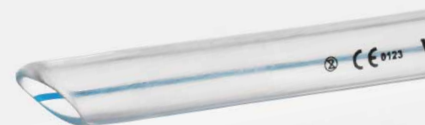

Order information

Gastric Tube Guide

for single use, Material: PVC, sterile

REF 49-40-018-1

Box 10

Specifications

O.D.	10.3 mm	Length	33.0 cm
I.D.	7.5 mm	Gastric Tube	max. 18 Fr

VBM Medizintechnik GmbH

Einsteinstrasse 1 | 72172 Sulz a. N. | Germany

Tel.: +49 7454 9596-0 | Fax: +49 7454 9596-33 | e-mail: info@vbm-medical.de | www.vbm-medical.de

